

Please type a plus sign (+) inside this box → ☐

Approved for use through 09/30/2000. OMB 0851-0032  
Patent and Trademark Office: U.S. DEPARTMENT OF COMMERCE

Under the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it displays a valid OMB control number.

# UTILITY PATENT APPLICATION TRANSMITTAL

(Only for new nonprovisional applications under 37 C.F.R. § 1.53(b))

Attorney Docket No. 3598-6  
First Inventor or Application Identifier Suk-song Oh  
Title Endodontic Gutta-percha With Working Length Marks  
Express Mail Label No.

**APPLICATION ELEMENTS**  
See MPEP chapter 600 concerning utility patent application contents.

**ADDRESS TO:** Assistant Commissioner for Patents  
Box Patent Application  
Washington, DC 20231

1. ☒ \* Fee Transmittal Form (e.g., PTO/SB/17)  
(Submit an original and a duplicate for fee processing)
2. ☒ Specification [Total Pages 12]  
(preferred arrangement set forth below)  
- Descriptive title of the invention  
- Cross References to Related Applications  
- Statement Regarding Fed sponsored R & D  
- Reference to Microfiche Appendix  
- Background of the invention  
- Brief Summary of the invention  
- Brief Description of the Drawings (if filed)  
- Detailed Description  
- Claim(s)  
- Abstract of the Disclosure
3. ☒ Drawing(s) (35 U.S.C. 113) [Total Sheets 4]  
4. Oath or Declaration [Total Pages 2]  
a. ☒ Newly executed (original or copy)  
b. ☐ Copy from a prior application (37 C.F.R. § 1.63(d))  
(for continuation/divisional with Box 17 completed)  
[Note Box 5 below]  
i. ☐ **DELETION OF INVENTOR(S)**  
Signed statement attached deleting  
inventor(s) named in the prior application,  
see 37 C.F.R. §§ 1.63(d)(2) and 1.33(b).
5. ☐ Incorporation By Reference (useable if Box 4b is checked)  
The entire disclosure of the prior application, from which a  
copy of the oath or declaration is supplied under Box 4b, is  
considered to be part of the disclosure of the accompanying  
application and is hereby incorporated by reference therein.

6. ☐ Microfiche Computer Program (Appendix)
7. Nucleotide and/or Amino Acid Sequence Submission  
(if applicable, all necessary)  
a. ☐ Computer Readable Copy  
b. ☐ Paper Copy (identical to computer copy)  
c. ☐ Statement verifying identity of above copies

## ACCOMPANYING APPLICATION PARTS

8. ☐ Assignment Papers (cover sheet & document(s))
9. ☐ 37 C.F.R. § 3.73(b) Statement (when there is an assignee) ☐ Power of Attorney
10. ☐ English Translation Document (if applicable)
11. ☒ Information Disclosure Statement (IDS)/PTO-1449 ☒ Copies of IDS Citations
12. ☐ Preliminary Amendment
13. ☒ Return Receipt Postcard (MPEP 503)  
(Should be specifically itemized)
14. ☒ \* Small Entity Statement(s) ☐ Statement filed in prior application,  
(PTO/SB/09-12) Status still proper and desired
15. ☒ Certified Copy of Priority Document(s)  
(if foreign priority is claimed)
16. ☐ Other: .....

\* NOTE FOR ITEMS 1 & 14: IN ORDER TO BE ENTITLED TO PAY SMALL ENTITY  
FEES, A SMALL ENTITY STATEMENT IS REQUIRED (37 C.F.R. § 1.27), EXCEPT  
IF ONE FILED IN A PRIOR APPLICATION IS RELIED UPON (37 C.F.R. § 1.28).

17. If a CONTINUING APPLICATION, check appropriate box, and supply the requisite information below and in a preliminary amendment:

☐ Continuation ☐ Divisional ☐ Continuation-in-part (CIP) of prior application No: \_\_\_\_\_  
Prior application information: Examiner \_\_\_\_\_ Group / Art Unit: \_\_\_\_\_

## 18. CORRESPONDENCE ADDRESS

☐ Customer Number or Bar Code Label (Insert Customer No. or Attach bar code label here) or ☒ Correspondence address below

Name M. Arthur Auslander  
Auslander & Thomas  
Address 505 Eighth Avenue  
City New York State NY Zip Code 10018  
Country U.S.A. Telephone 212-594-6900 Fax 212-244-0028

Name (Print/Type) M. Arthur Auslander Registration No. (Attorney/Agent) 18,376  
Signature [Signature] Date April 16, 1999

Burden Hour Statement: This form is estimated to take 0.2 hours to complete. Time will vary depending upon the needs of the individual case. Any  
comments on the amount of time you are required to complete this form should be sent to the Chief Information Officer, Patent and Trademark Office,  
Washington, DC 20231. DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADDRESS. SEND TO: Assistant Commissioner for Patents,  
Box Patent Application, Washington, DC 20231.

Under the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it displays a valid OMB control number.

**STATEMENT CLAIMING SMALL ENTITY STATUS  
(37 CFR 1.9(f) & 1.27(b))--INDEPENDENT INVENTOR**Docket Number (Optional)  
3598-6Applicant, Patentee, or Identifier: Suk-song Oh

Application or Patent No.: \_\_\_\_\_

Filed or Issued: \_\_\_\_\_

Title: Endodontic Gutta-percha Point With Working Length Marks

As a below named inventor, I hereby state that I qualify as an independent inventor as defined in 37 CFR 1.9(c) for purposes of paying reduced fees to the Patent and Trademark Office described in:

- ☒ the specification filed herewith with title as listed above.
- ☐ the application identified above.
- ☐ the patent identified above.

I have not assigned, granted, conveyed, or licensed, and am under no obligation under contract or law to assign, grant, convey, or license, any rights in the invention to any person who would not qualify as an independent inventor under 37 CFR 1.9(c) if that person had made the invention, or to any concern which would not qualify as a small business concern under 37 CFR 1.9(d) or a nonprofit organization under 37 CFR 1.9(e).

Each person, concern, or organization to which I have assigned, granted, conveyed, or licensed or am under an obligation under contract or law to assign, grant, convey, or license any rights in the invention is listed below:

- ☐ No such person, concern, or organization exists.
- ☒ Each such person, concern, or organization is listed below.


Suk-song Oh  
426-11 Shinsu-dong, Mapo-gu  
Seoul, Korea

Separate statements are required from each named person, concern, or organization having rights to the invention stating their status as small entities. (37 CFR 1.27)

I acknowledge the duty to file, in this application or patent, notification of any change in status resulting in loss of entitlement to small entity status prior to paying, or at the time of paying, the earliest of the issue fee or any maintenance fee due after the date on which status as a small entity is no longer appropriate. (37 CFR 1.28(b))

Suk-song, OH

NAME OF INVENTOR

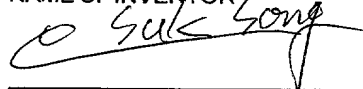

Signature of inventor

Date

1999. 4. 6.

NAME OF INVENTOR

Signature of inventor

Date

NAME OF INVENTOR

Signature of inventor

Date

# ENDODONTIC GUTTA-PERCHA POINT WITH WORKING LENGTH MARKS

## BACKGROUND OF THE INVENTION

### Field of the Invention

The present invention relates, in general, to gutta-percha points used as plugging materials for root canals in endodontic treatments and, more particularly, to an endodontic gutta-percha point having working length marks at its upper portion.

### Description of the Prior Art

As well known to those skilled in the art, gutta percha is a material, which is extracted from palaguium or payena generic plants of sapotaceae, the native plants of the Malay Archipelago. The gutta percha has a semisolid phase at an ordinary temperature, but is changed in the phase when it is pressurized or heated. When the gutta percha is heated to a temperature higher than 60°C, it becomes a plastic.

Gutta-percha points, used as plugging materials for a root canal in an endodontic treatment, may be different in their concentrations of components according to manufacturers of the points. However, most of gutta-percha points typically consist of 66% of zinc oxide, 20% of gutta percha, 11% of heavy metal sulphate (a radio-opacity material), and 3% of wax or resin used as a plasticizer.

Such gutta-percha points are problematic in that it is somewhat difficult to insert the points into narrow root canals and it is necessary to use a sealer in order to seal the root canals after inserting the points into the canals.

Another problem, experienced in typical gutta-percha points, resides in that they may be easily broken when a period of time has passed after manufacture.

However, the gutta-percha points preferably meet most of the needs typically desired from plugging materials for root canals in endodontic treatments.

Since the gutta-percha points are plastic materials, they are preferably accepted by an access cavity for the root canal and are easy to treat.

In addition, it is somewhat easy to remove the gutta-percha points from root canals in the case of a root canal retreatment or a dentodic operation. The gutta-percha points also have a low toxic effect to an apical tissue.

Gutta-percha points, which are marketed now, are generally classified into two types: standardized points and conventional points.

The standardized gutta-percha points are sized and tapered in the same scale as that of typical endodontic instruments, thus being more generally used by endodontists.

Such gutta percha has been used as a material of endodontic plugging points for 100 years or more. However, there has not been proposed another material which is more acceptable or preferable than the gutta percha as a material of the endodontic plugging points. Therefore, the gutta-percha points are widely and preferably used as root canal plugging points in endodontic treatments.

Figs. 3a to 3g schematically show an endodontic treatment process of treating a decayed root canal using the conventional gutta-percha points.

A root canal is typically infected as follows. As shown

in Fig. 3a, a tooth "T" is initially decayed at the enamel 31. The decaying of the tooth is, thereafter, propagated to the dentine 32, the pulp cavity 33 and the root canal 34 in order. The pulp tissue in the tooth T is also infected in accordance with the propagation of decaying. When the pulp tissue is decayed, a patient feels pain.

A root canal treatment for such a decayed tooth T is carried out as follows. In order to treat the decayed tooth T, it is primarily necessary to prepare an access cavity 35 in the tooth T as shown in Fig. 3b.

In order to form such an access cavity 35, the decayed portions of both the enamel 31 and the dentine 32 are removed from the tooth T using an endodontic instrument, such as a file or a reamer. In such a case, the decayed pulp tissue has to be also removed. The access cavity 35 is enlarged by smoothing the interior surface of the cavity 35 using the endodontic instrument. The enlarged and smoothed cavity 35 is finally sterilized and dried using paper points.

After the access cavity 35 is primarily formed in the tooth T, it is necessary to measure a working length WL from the top of the tooth T to the apical foramen 36 as shown in Fig. 3c.

In such a case, the measurement of the working length WL in the decayed tooth T is for limiting the access area for endodontic instruments, medicines and plugging materials so as to prevent the apical tissue from any damage during both an access cavity preparation process and a root canal plugging process.

Of course, there may be different working lengths WL in accordance with patients, teeth T, and positions of apical

foramens 36. However, such working lengths WL typically range from about 18 mm to about 24 mm.

The measurement of the working length WL is performed using a file 37 with a silicon rubber stopper 38. That is, in order to measure the working length WL, the file 37 is primarily inserted into the access cavity 35 until the apex of the file 37 reaches a position around the apical foramen 36 without piercing the apical tissue or giving pain to the tissue. Thereafter, an X-ray inspection is carried out so as to determine whether the apex of the file 37 is appropriately accepted by the apical foramen 36 or not. When the X-ray inspection shows that the apex of the file 37 is appropriately accepted by the apical foramen 36, the stopper 38 is adjusted along the file 37 until the stopper 38 meets the top of the tooth T. In such a case, the length from the apex of the file 37 to the stopper 38 is the working length WL.

The working length WL on the file 37 is measured using an endodontic ruler which shows the length WL as a numerical value. Thereafter, gutta-percha points, including master and accessory points, are measured using the endodontic ruler so as to determine a length, which is equal to the calculated working length WL, on each point. In such a case, the length measurement of the gutta-percha points is started from the sharpened apex of each point. After measuring a working length WL of each gutta-percha point using the endodontic ruler, the end of the length WL opposite to the apex is marked on each point prior to starting a root canal plugging process using the points. The root canal plugging process is shown in Fig. 3d.

In the root canal plugging process, a sealer is primarily

applied to the access cavity 35. Thereafter, the master point of the gutta-percha points is inserted into the tooth T until the mark of the working length WL on the master point is aligned with the top of the tooth T. The insertion of the master point is followed by an insertion of accessory points. That is, after the insertion of the master point, a plurality of accessory points of the gutta-percha points are inserted into the access cavity 35 at positions around the master point using a spreader. The insertion of the accessory point is also carried out through a lateral condensation method. In the root canal plugging process, it is also necessary to prevent any gap from being left in the access cavity 35 filled with the points.

In the above plugging process, it is preferable to make the apex of each gutta-percha point precisely aligned with the apical foramen 36 of the access cavity 35.

When the apex of each gutta-percha point fails to reach the apical foramen 36, a dead space is formed at an area between the apex of the point and the apical foramen 36. Such a dead space may allow a propagation of bacteria at that space and this forces the patient to suffer from the root canal retreatment. On the other hand, when the sharpened apex of each gutta-percha point exceeds the apical foramen 36 and pierces into the apical tissue, the apex of the point pressurizes the apical tissue and gives the patient pain.

Therefore, it is necessary to insert the gutta-percha points, particularly, the master point, into the access cavity 35 with the apex of each point being appropriately accepted by the apical foramen 36. In order to accomplish the above object, it is very important to precisely mark the working

length WL on each point, particularly, the master point, and to precisely insert the points into the access cavity 35 by the working length WL.

After inserting the gutta-percha points into the access cavity 35, the upper portions of the points are removed using a plugger with the lower portions of the points being left in the enlarged root canal as shown in Fig. 3e.

Fig. 3f shows an amalgam 39 filling in the upper portion of the access cavity 35 with the gutta-percha points plugging the enlarged root canal.

As shown in Fig. 3g, a gold crown 40 is, thereafter, exteriorly set on the top portion of the tooth T filled with both the gutta-percha points and the amalgam 39 in the access cavity 35. An endodontic treatment process is thus accomplished.

However, the above-mentioned process gives inconvenience to both endodontists and patients due to a structural defect of the typical gutta-percha points. That is, in the above process, it is necessary for an endodontist to measure a working length WL using a file 37. The working length WL on the file 37 is measured using an endodontic ruler which shows the length WL as a numerical value. Thereafter, gutta-percha points, including master and accessory points, are measured using the endodontic ruler so as to determine a length, which is equal to the calculated working length WL, on each point prior to marking the upper end of the length WL opposite to the apex of each point. However, typical gutta-percha points are free from any working length marks on its surface as shown in Fig. 1. Therefore, it is necessary for endodontists to mark the working lengths WL on the gutta-percha points



during root canal treatments. In order to mark the working lengths WL on the points 20, each point 20 may be bent at a position corresponding to the upper end of the length WL, thus forming a bent portion 21 on the point 20 as shown in Fig. 2a. Alternatively, the working length WL may be marked on each point 20 by pressing the point 20 using an appropriate marking instrument, such as a pincette, at a position corresponding to the upper end of the length WL, thus forming a pressed mark 22 on the point 20. Of course, there may be other methods of forming the mark of the working length WL on each of the gutta-percha points 20.

Therefore, the typical gutta-percha points 20 are problematic in that they force an endodontist to repeatedly and manually measure a working length WL and to form the bent portions 21 or the pressed marks 22 on the points 20 during a root canal treatment. Such gutta-percha points 20 thus consume time and this is inconvenient to both the endodontists and patients since it is necessary to precisely perform a root canal treatment within a somewhat short period of time with the mouth of a patient being continuously and fully opened.

Furthermore, there may occur an undesired error while measuring a working length WL with an endodontic ruler or forming a bent portion 21 or a pressed mark 22 on each point 20 so as to mark the working length WL on each point 20 during a root canal treatment. Such an error makes the apexes of the points 20 fail to be precisely aligned with an apical foramen 36. This regrettably results in an operational error during a root canal treatment.

#### SUMMARY OF THE INVENTION

Accordingly, the present invention has been made keeping in mind the above problems occurring in the prior art, and an object of the present invention is to provide an endodontic gutta-percha point having working length marks at its upper portion.

In order to accomplish the above object, the present invention provides an endodontic gutta-percha point, which has a plurality of engraved or embossed marks as working length marks on its upper portion, thus allowing endodontists to be free from a repeated marking of working lengths on the gutta-percha points while plugging an enlarged root canal with the points during an endodontic treatment. This allows the root canal plugging process to be accomplished within a desired short period of time. The working length marks of the gutta-percha points also allow an endodontist to precisely insert the points into an access cavity of a decayed tooth by a desired depth without failure.

#### BRIEF DESCRIPTION OF THE DRAWINGS

The above and other objects, features and other advantages of the present invention will be more clearly understood from the following detailed description taken in conjunction with the accompanying drawings, in which:

Fig. 1 is a front view of a conventional gutta-percha point used as a root canal plugging material in an endodontic treatment;

Figs. 2a and 2b are views, respectively showing typical methods of marking a working length on a conventional gutta-percha point;

Figs. 3a to 3g are views, schematically showing an endodontic treatment process of treating a decayed root canal using such conventional gutta-percha points;

Fig. 4 is a front view of a gutta-percha point with a plurality of working length marks in accordance with the primary embodiment of the present invention;

Fig. 5 is a view showing the gutta-percha point of Fig. 4 which is appropriately inserted into an access cavity of a decayed tooth during an endodontic treatment process; and

Fig. 6 is a front view of a gutta-percha point with a plurality of working length marks in accordance with the second embodiment of the present invention.

#### DESCRIPTION OF THE PREFERRED EMBODIMENTS

Fig. 4 is a front view of a gutta-percha point with a plurality of working length marks in accordance with the primary embodiment of this invention. As shown in the drawing, the gutta-percha point 1 according to the invention has a plurality of working length marks 2 along its upper portion, with the marks 2 being spaced apart from each other. In the primary embodiment, the marks 2 are formed on the point 1 through an embossing process.

The working length marks 2 are formed on the upper portion of the point 1 within a marking range, the range starting at a position spaced apart from the sharpened apex of the point 1 by 18 mm and ending at a position spaced apart from said apex by 24 mm. The marks 2 are also spaced out at different intervals of 1 mm and 2 mm.

In order to allow the marks 2 to be somewhat prominently

distinguished by users or endodontists, it is preferable to paint the circumferential surface of each embossed mark 2 with a distinguishable color.

The above gutta-percha point 1 is used as follows. Fig. 5 shows the point 1 appropriately inserted into an access cavity of a decayed tooth during an endodontic treatment process. When the point 1 is used for plugging an enlarged root canal of a decayed tooth T having a working length WL of 18 mm, the point 1 is inserted into an access cavity 35 of the tooth T until an 18 mm mark 2, corresponding to the 18 mm working length, is precisely aligned with the top of the tooth T. When the 18 mm mark 2 is precisely aligned with the top of the tooth T, the apex of the point 1 is precisely aligned with the apical foramen 36, thus being free from forming a dead space in the root canal or passing through the apical foramen 36 so as to pierce the apical tissue.

Fig. 6 shows a gutta-percha point with a plurality of working length marks in accordance with the second embodiment of this invention. In the second embodiment, the gutta-percha point 11 has a plurality of working length marks 12 along its upper portion in the same manner as that described for the primary embodiment. However, the marks 12 of the second embodiment are formed on the point 11 through an engraving process different from the primary embodiment.

In the same manner as that described for the primary embodiment, the working length marks 12 are formed on the upper portion of the point 11 within a marking range. The marking range starts at a position spaced apart from the sharpened apex of the point 11 by 18 mm and ends at a position spaced apart from said apex by 24 mm. The marks 12 are also

spaced out at different intervals of 1 mm and 2 mm.

In order to allow the marks 12 to be somewhat prominently distinguished by users or endodontists, the circumferential surface of each embossed mark 2 is preferably painted with a distinguishable color.

The above gutta-percha point 11 is used as follows. When the point 11 is used for plugging an enlarged root canal of a decayed tooth T having a working length WL of 18 mm, the point 11 is inserted into an access cavity 35 of the tooth T until an 18 mm mark 12, corresponding to the 18 mm working length, is precisely aligned with the top of the tooth T. When the 18 mm mark 12 is precisely aligned with the top of the tooth T, the apex of the point 1 is precisely aligned with the apical foramen 36, thus being free from forming a dead space in the root canal or passing through the apical foramen 36 so as to pierce the apical tissue.

As described above, the present invention provides a gutta-percha point used as a root canal plugging material in an endodontic treatment. The gutta-percha point of this invention has a plurality of engraved or embossed marks as working length marks on its upper portion. The point thus allows endodontists to be free from a repeated marking of working lengths on the gutta-percha points while plugging an enlarged root canal with the points during an endodontic treatment. This allows the root canal plugging process to be accomplished within a desired short period of time. The working length marks of the gutta-percha points also allow an endodontist to precisely insert the points into an access cavity of a decayed tooth by a desired depth without failure. Therefore, the gutta-percha point is free from forming a dead

space in the root canal. or passing through an apical foramen  
so as to pierce an apical tissue. The gutta-percha points of  
this invention thus allow endodontists to appropriately  
accomplish desired root canal treatments within a short period  
of time.

Although the preferred embodiments of the present  
invention have been disclosed for illustrative purposes, those  
skilled in the art will appreciate that various modifications,  
additions and substitutions are possible, without departing  
from the scope and spirit of the invention as disclosed in the  
accompanying claims.

WHAT IS CLAIMED IS:

1. A gutta-percha point used as a root canal plugging material in an endodontic treatment, comprising:

5 a plurality of working length marks along an upper portion of the point, said marks being spaced apart from each other.

10 2. The gutta-percha point according to claim 1, wherein said working length marks are formed through an embossing process.

15 3. The gutta-percha point according to claim 1, wherein said working length marks are formed through an engraving process.

20 4. The gutta-percha point according to claim 2, wherein each of the embossed working length marks is painted with a distinguishable color along a circumferential surface thereof.

5. The gutta-percha point according to claim 3, wherein each of the engraved working length marks is painted with a distinguishable color along a circumferential surface thereof.

25 6. The gutta-percha point according to any one of claims 2 to 5, wherein said working length marks are formed on the upper portion of said point within a marking range and are spaced out at intervals of 1 mm and 2 mm, said marking range starting at a position spaced apart from an apex of the point  
30 by 18 mm and ending at a position spaced apart from said apex

by 24 mm.




## ABSTRACT

A gutta-percha point, used as a root canal plugging material in an endodontic treatment, is disclosed. The gutta-percha point has a plurality of engraved or embossed marks as working length marks on its upper portion. Each of the engraved or embossed marks is preferably painted with a distinguishable color along a circumferential surface. The above marks are formed on the point within a marking range while being spaced out at intervals of 1 mm and 2 mm. The marking range starts at a position, spaced apart from the apex of the point by 18 mm, and ends at a position spaced apart from the apex by 24 mm. The gutta-percha point thus allows an endodontist to be free from a repeated marking of working lengths on such points during an endodontic treatment, thus allowing the endodontic treatment to be accomplished within a desired short period of time. Due to the working length marks, the gutta-percha point is precisely inserted into an access cavity of a decayed tooth by a desired depth, thus being free from forming a dead space in the root canal or piercing the apical tissue.


**FIG. 1**

Prior Art

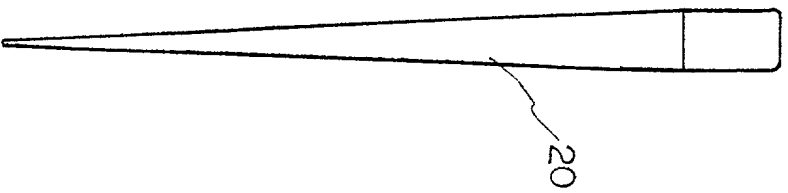


**FIG. 2a**

Prior Art

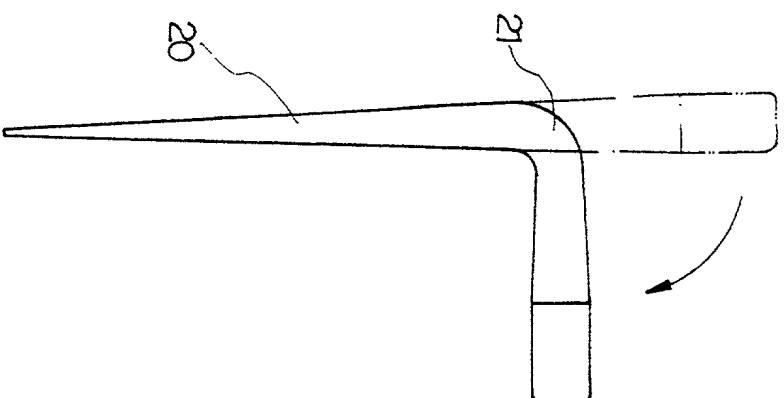


**FIG. 2b**


Prior Art

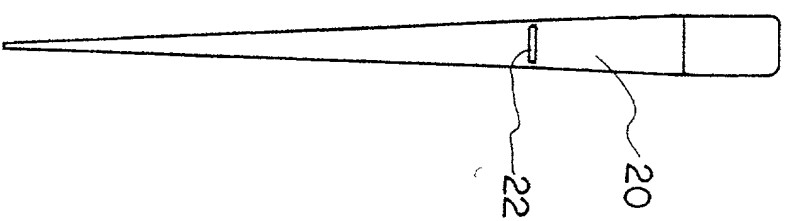


**FIG. 3a**

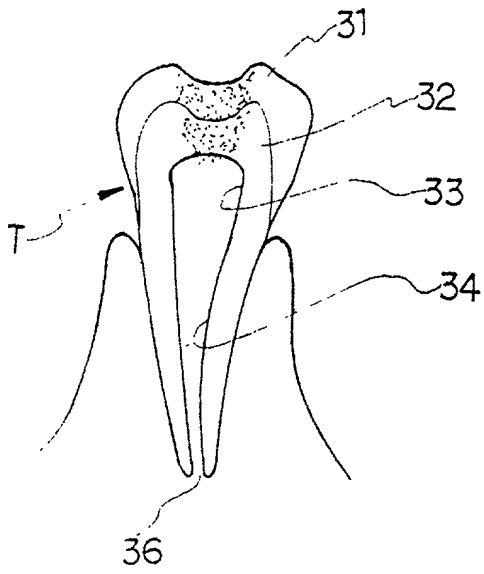


**FIG. 3b**

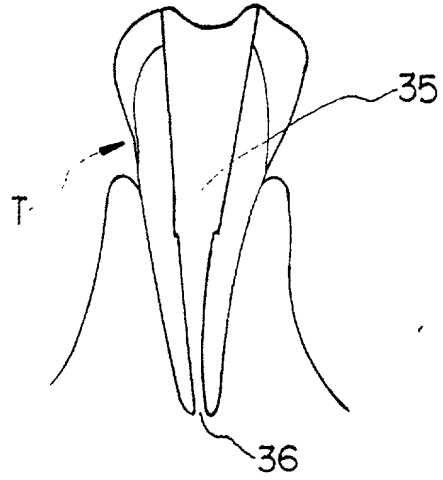


**FIG. 3c**

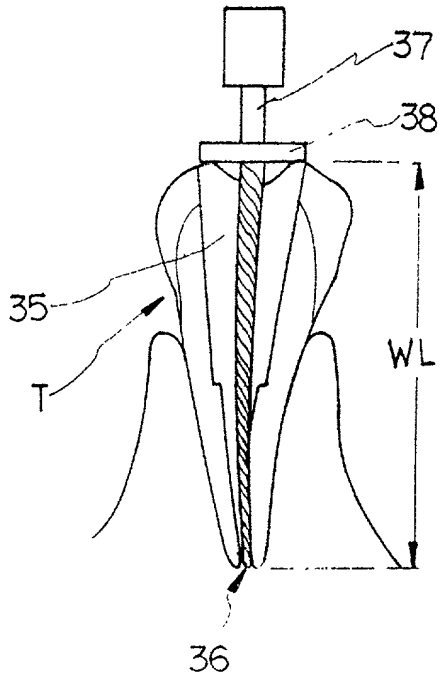


**FIG. 3d**

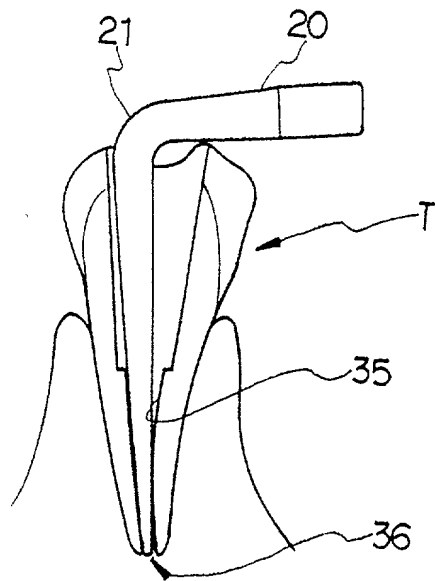

**FIG. 3e**

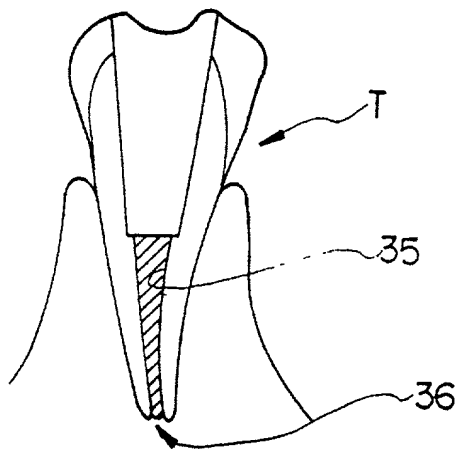

**FIG. 3f**

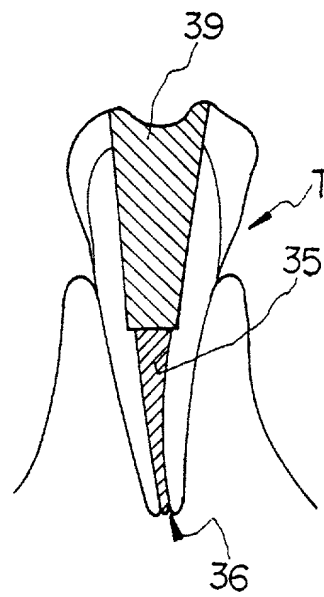

**FIG. 3g**

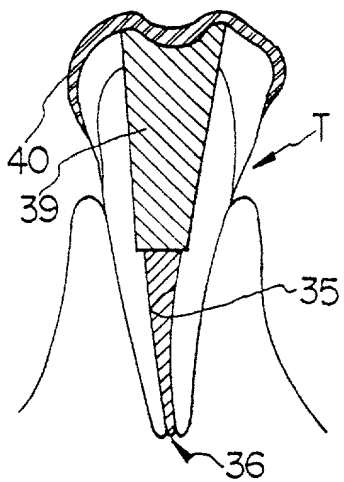


FIG.4

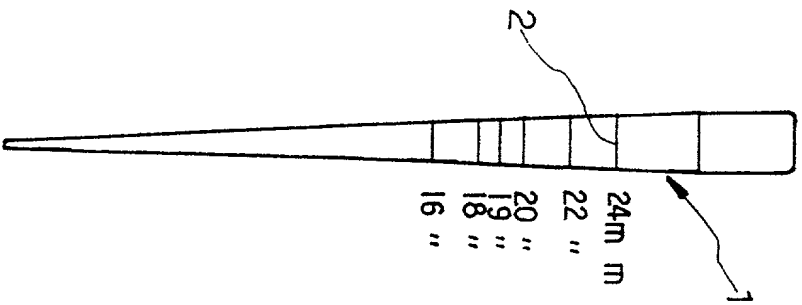


FIG.5


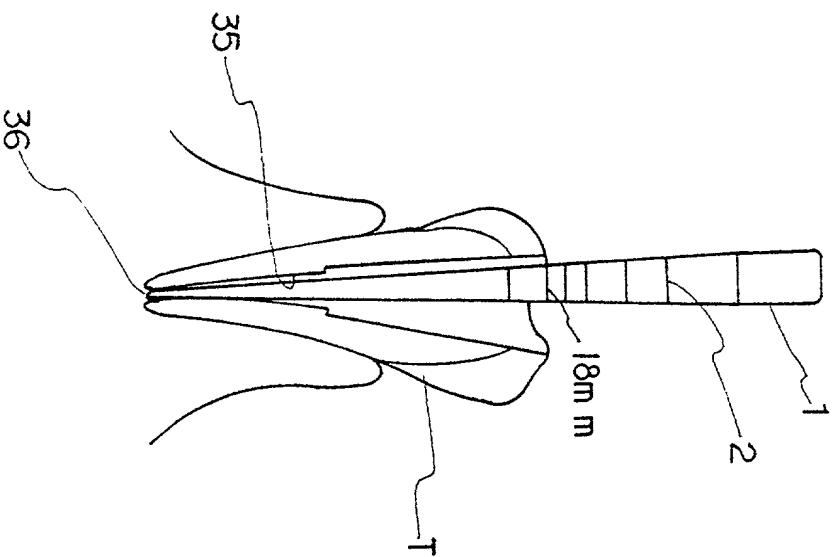


FIG.6

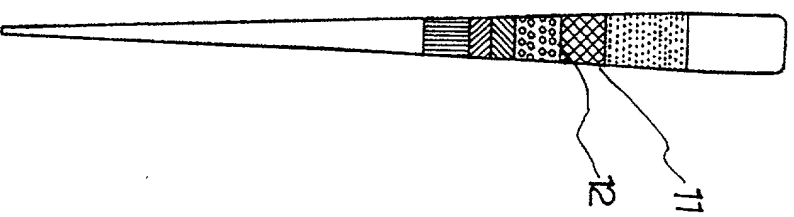

Under the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it contains a valid OMB control number.

<b>DECLARATION FOR UTILITY OR DESIGN PATENT APPLICATION (37 CFR 1.63)</b>	<b>Attorney Docket Number</b>		3598-6
	<b>First Named Inventor</b>		Suk-song Oh
	<b>COMPLETE IF KNOWN</b>		
	<b>Application Number</b>	/	
	<b>Filing Date</b>		
	<b>Group Art Unit</b>		
<b>Examiner Name</b>			

☒ Declaration Submitted with Initial Filing
 OR
 ☐ Declaration Submitted after Initial Filing (surcharge (37 CFR 1.16 (e)) required)

**As a below named inventor, I hereby declare that:**

My residence, post office address, and citizenship are as stated below next to my name.

I believe I am the original, first and sole inventor (if only one name is listed below) or an original, first and joint inventor (if plural names are listed below) of the subject matter which is claimed and for which a patent is sought on the invention entitled:

## Endodontic Gutta-percha Point With Working Length Marks

the specification of which

*(Title of the Invention)*

☒ is attached hereto  
OR

☐ was filed on (MM/DD/YYYY) [redacted] as United States Application Number or PCT International

Application Number  and was amended on (MM/DD/YYYY)  (if applicable).

I hereby state that I have reviewed and understand the contents of the above identified specification, including the claims, as amended by any amendment specifically referred to above.

I acknowledge the duty to disclose information which is material to patentability as defined in 37 CFR 1.56.

I hereby claim foreign priority benefits under 35 U.S.C. 119(a)-(d) or 365(b) of any foreign application(s) for patent or inventor's certificate, or 365(a) of any PCT international application which designated at least one country other than the United States of America, listed below and have also identified below, by checking the box, any foreign application for patent or inventor's certificate, or of any PCT international application having a filing date before that of the application on which priority is claimed.

Prior Foreign Application Number(s)	Country	Foreign Filing Date (MM/DD/YYYY)	Priority Not Claimed	Certified Copy Attached?	
				YES	NO
1998-37719	Korea	September 12, 1998	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>
			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

☐ Additional foreign application numbers are listed on a supplemental priority data sheet PTO/SB/02B attached hereto:

I hereby claim the benefit under 35 U.S.C. 119(e) of any United States provisional application(s) listed below.

Application Number(s)	Filing Date (MM/DD/YYYY)	
		<input type="checkbox"/> Additional provisional application numbers are listed on a supplemental priority data sheet PTO/SB/02B attached hereto.

[Page 1 of 2]

**Burden Hour Statement:** This form is estimated to take 0.4 hours to complete. Time will vary depending upon the needs of the individual case. Any comments on the amount of time you are required to complete this form should be sent to the Chief Information Officer, Patent and Trademark Office, Washington, DC 20231. **DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADDRESS. SEND TO: Assistant Commissioner for Patents, Washington, DC 20231.**

Please type a plus sign (+) inside this box → ☐

PTO/SB/01 (12-97)

Approved for use through 9/30/00. OMB 0651-0032

Patent and Trademark Office; U.S. DEPARTMENT OF COMMERCE  
Under the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it contains a valid OMB control number.

## DECLARATION — Utility or Design Patent Application

I hereby claim the benefit under 35 U.S.C. 120 of any United States application(s), or 365(c) of any PCT international application designating the United States of America, listed below and, insofar as the subject matter of each of the claims of this application is not disclosed in the prior United States or PCT International application in the manner provided by the first paragraph of 35 U.S.C. 112, I acknowledge the duty to disclose information which is material to patentability as defined in 37 CFR 1.56 which became available between the filing date of the prior application and the national or PCT international filing date of this application.

U.S. Parent Application or PCT Parent Number	Parent Filing Date (MM/DD/YYYY)	Parent Patent Number (if applicable)

☐ Additional U.S. or PCT international application numbers are listed on a supplemental priority data sheet PTO/SB/02B attached hereto.

As a named inventor, I hereby appoint the following registered practitioner(s) to prosecute this application and to transact all business in the Patent and Trademark Office connected therewith:

☐ Customer Number  OR  
☐ Registered practitioner(s) name/registration number listed below

Place Customer Number Bar Code Label here

Name	Registration Number	Name	Registration Number
M. Arthur Auslander	18,376		


☐ Additional registered practitioner(s) named on supplemental Registered Practitioner Information sheet PTO/SB/02C attached hereto.

Direct all correspondence to: ☐ Customer Number  OR ☒ Correspondence address below

Name	M. Arthur Auslander				
Address	Auslander & Thomas				
Address	505 Eighth Avenue				
City	New York	State	NY	ZIP	10018
Country	U.S.A.	Telephone	212-594-6900	Fax	212-244-0028

I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under 18 U.S.C. 1001 and that such willful false statements may jeopardize the validity of the application or any patent issued thereon.

Name of Sole or First Inventor: ☐ A petition has been filed for this unsigned inventor

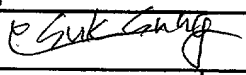Given Name (first and middle (if any))		Family Name or Surname			
Suk-song		Oh			
Inventor's Signature				Date	1999. 4. 6.
Residence: City	Seoul	State		Country	Korea
				Citizenship	Rep. of Korea
Post Office Address	426-11 Shinsu-dong, Mapo-gu				
Post Office Address					
City	Seoul	State		ZIP	
				Country	Korea

☐ Additional inventors are being named on the \_\_\_\_\_ supplemental Additional Inventor(s) sheet(s) PTO/SB/02A attached hereto